


KODAK Precise[®] PB Progressive Lens Dispensing Aid



Kodak
Precise PB Lens

Kodak
Precise Short PB Lens

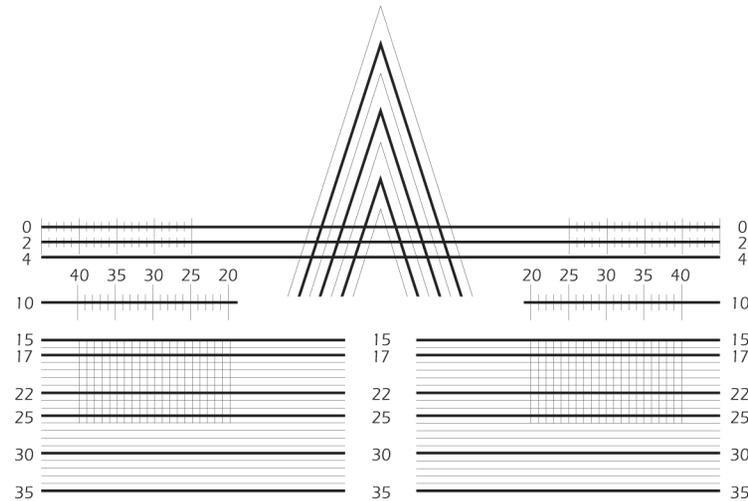

KODAK Precise PB Progressive Lens

1. Select and adjust the frame. Adjust the frame for comfort and accuracy before taking measurements. Adjustable nose pads are recommended. Set the vertex distance 12-14mm. Set the pantoscopic angle to 10-12°. Frame should have a slight face form.

2. PD and Fit Height. Measure monocular fitting heights by marking each demo lens at the pupil centers with a felt tip pen. Measure monocular PDs using a pupilometer or by using the fitting height marks. To translate lens markings into measurements, use the scale shown above. Place the frame on the center of the triangle, making sure the marks on the lens are on the zero (0) line. Using the scale, record the monocular PDs and monocular vertical seg heights.

(Continue to remaining steps 3 through 5 on the reverse side of this page.)


Signet Armorlite, Inc. 5803 Newton Drive, Suite A
Carlsbad, CA 92008
800-759-0075
www.signetarmorlite.com

To access marketing materials, log on to www.SALitOnline.com

See the *Colors* of Life[®]

KODAK Precise PB Progressive Lens Dispensing Aid

KODAK Precise PB Progressive Lens



3. Frame verification. Line up the pupillary mark on the frame with the cross on the lens chart to the right. **The Blue Optical Zone** is the area required for good distance and near vision. Verify that it is within the eyewire and that the eyewire is within the cut-out diameter. If a small portion of the near or distance zones are outside the frame, the patient may still be satisfied with the frame, as long as the optical compromises are explained.

4. Dispensing. Confirm the monocular PDs and fitting heights. Verify the prescription. Confirm the fit on the patient by verifying that the fitting cross is properly positioned over the pupils. Adjust the frame as necessary.

5. Teach proper viewing. Demonstrate the different viewing areas, and appropriate head and eye movement.

KODAK Precise PB Lenses



3.0mm Decentered

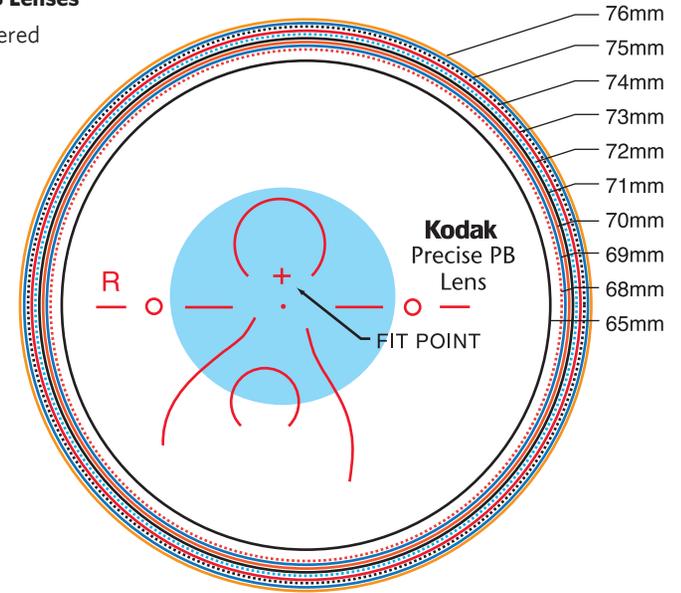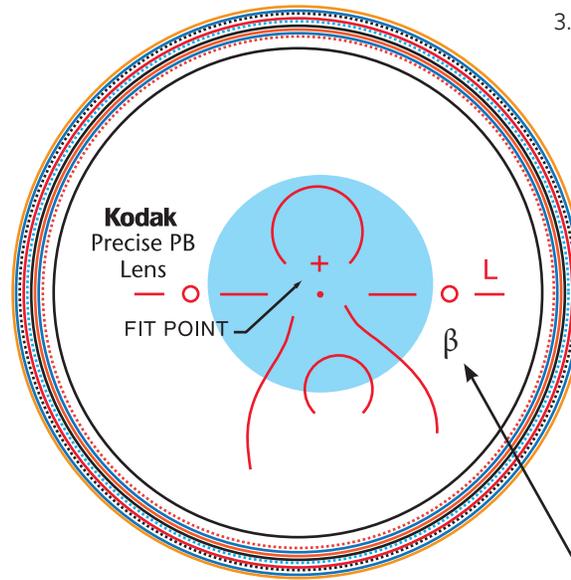

KODAK Precise PB Lens identifying mark

KODAK Precise Short PB Lenses

3.0mm Decentered

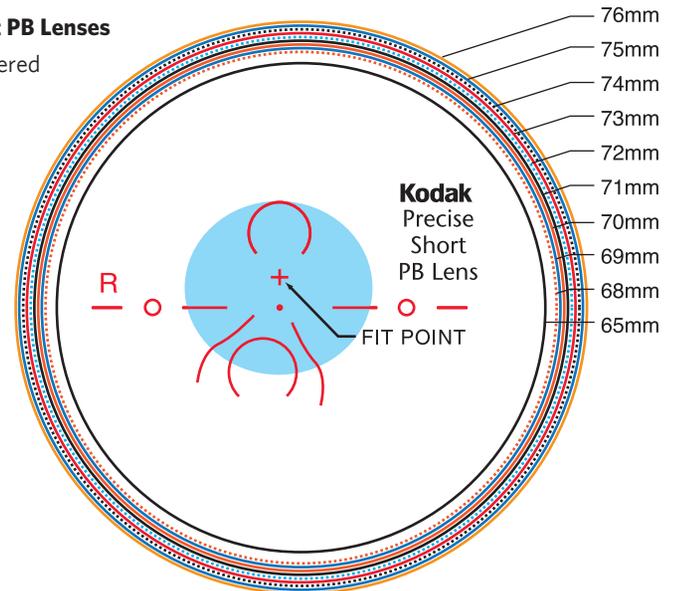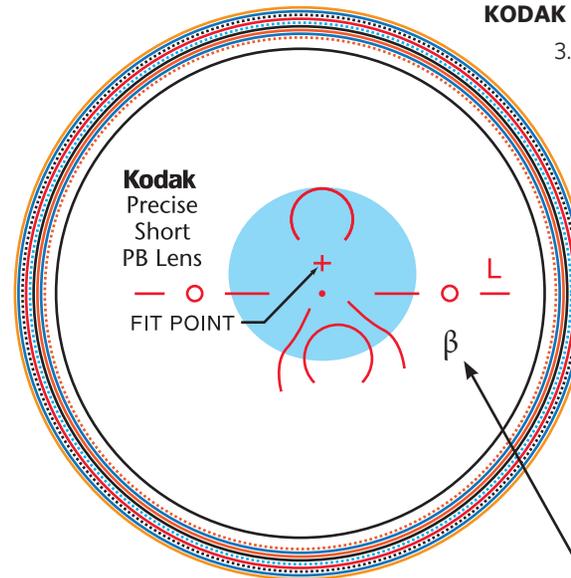

KODAK Precise PB Lens identifying mark



Exclusively available to
Independent Eyecare Professionals
Learn more at www.kodaklens.com/pro